

CJ-Optik GmbH & Co. KG

Willeckstrasse 1 35614 Asslar-Werdorf Germany

Phone: +49 (0) 64 43 - 81 81 33 0 +49 (0) 64 43 - 81 81 33 9 Fax:

info@cj-optik.de www.cj-optik.de





TajhizTeb Asar Co. Exclusive agent of CJ-Optik in Iran.

Unit 35, No. 1831 North Kargar Blvd. Tehran, Iran.

Phone: + 98 21 42 86 8000 + 98 21 42 86 8102 Fax: info@tajhiztebasar.com

www.tajhiztebasar.com













Advanced & Advanced SensorUnit



German Design Award 2017 for Flexion







Statement of the jury:

"A high-performance high-tech product with a very progressive and contemporary design."

Certificates / Awards











Advanced





The Flexion Family







Basic





Flexion Advanced / Flexion Advanced SensorUnit

Wetzlar - City of Optics













Many world-famous companies from the optics sector have settled in Wetzlar since 1850. Today, over 70 optics companies are located in the city and the surrounding area. CJ-Optik takes advantage of the region's outstanding expertise by cooperating with the "Wetzlar Network" and local suppliers.

Timeline

2007

2008

2013

2014

2016

Foundation of CJConsultation in Wetzlar



Developing of various OEM products



Foundation of CJ-Optik in Wetzlar



Introduction of VarioFocus 200-350 mm objective which fits for all microscopes



Relocation to Asslar. New factory building (1000 sqm)



2016

2017

2019

2021

2022

ISO 9001:2008 certification



German Design Award Special for the Flexion dental microscope



Production start of Flexion twin microscope



ISO 13485:2016 certification for medical devices



To date, over 4000 Flexion microscopes have been sold





Technical Data

	Flexion Advanced / Flexion Advanced SensorUnit			
Magnification System	5-step apochromatic magnification changer / 0.4 - 2.5			
Eyepieces	Widefield eyepieces WF 12.5 x Widefield eyepieces WF 12.5 x / Reticle *			
Tube System	$F = 170 \text{ mm apochromatic} / 180^{\circ} \text{ tiltable tube}$			
Stereo Base	24 mm for enhanced 3D vision			
Focus System	VarioFocus ² 200 - 350 mm / Plan-APO / glass protection * VarioFocus ³ 210 - 470 mm / Plan-APO / glass protection * VarioFocus with hydrophobic glass protection *			
Illumination System	Fanless high power LED 5.400 K / 65.000 h life time 2° with stepless spot diaphragm system Up to 100.000 Lux			
Filter System	Swing in green filter Swing in orange filter			
Movement	MonoGlobe, 355° Axis control			
Handles	Ergonomic handles, 360° adjustable			
Asepsis	1 x Asepsis set for the handles and knobs (sterilizable)			
Documentation	Integrated beam splitter 50:50 4K-Imaging-Port / FullFrame * HD-Imaging-Port / APS-C * Phone-Imaging-Port / Smartphones * SensorUnit controls for photo / video (only with Flexion Advanced SensorUnit) MicVision App (iOS & Android) * MicroVis Software *			
Power	Power adapter for all Sony Alpha cameras, smartphones			
Stands	Mobile stand *, Wall mount *, Ceiling mount *, Floor mount *, Dental units *			

* optional






Flexion Advanced / Flexion Advanced SensorUnit

Made in Germany



Our in-house production is located in Asslar in the immediate vicinity of the "City of Optics" Wetzlar. The production is according to CE, ISO 13485 and ISO 9001:2015 thus guaranteeing the highest quality and reliability. We take advantage of the region's out-standing expertise by cooperating with the "Wetzlar Network" and local suppliers.



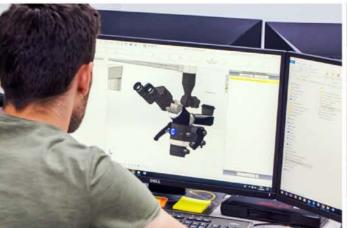








Accessories



Binoculars f 200

Compared to the normal focal length of 170 mm, the new f 200 binoculars lead to an increased magnification of 15 %.

Compatible with all CJ-Optik Flexion Microscopes.

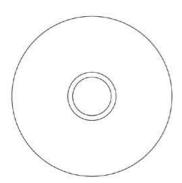

Eyepiece with Reticule

The reticule eyepiece allows you to center the positon of the image. This will help your assistance and for photo and video purposes.

Compatible with all CJ-Optik Flexion Microscopes.



Safety First: Laser Filter

We put safety first when you are working with a laser tool. Our laser filter with a range of 705 - 1290 nm ensures that your eyes are safe at any time. Easy assembly: All you have to do is remove the binocular head and then insert the filter. When the filter is not needed, simply pull it out a little.

Compatible with all CJ-Optik Flexion Microscopes.



Splash Guard Protection

CJ-Optik offers a splash guard to protect dentists and patients during treatment with the Flexion microscopes. With a simple assembly without tools, the plexiglass plates for the eyepiece and for the lens are attached. The set is available in two sizes.

Compatible with all CJ-Optik Flexion Microscopes.



Flexion Advanced

Since 2007 CJ-Optik has become one of the highperformance companies in the field of dental microscopy. With distributors in 75 countries worldwide we go for efficient solutions in the interests of dentists and patients.

Customer satisfaction is the cornerstone of our success. That is why we do everything we can to maintain customer confidence.

Top-quality products help us gain customer satisfaction. We can rely on many years of experience and technical progress. Our trained team ensures that our customers are looked after optimally at all times.


All over the world dentists and dental clinics share their practical experience with us. This led to the development of the new Flexion Advanced & Flexion Advanced SensorUnit microscopes.



Advanced Microscopy

Ergotube 30°

With binocular inclination 0 - 200°

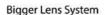

The larger the lenses and the better their coatings, the brighter and sharper the image

Large Exit Pupil

For wearers of glasses

IPD Adjuster

Single-handed, wide adjustment range from 55 - 75 mm

SensorUnit

Control photo or video recording easily with a hand gesture

Superior Light Transmission

Throughput for optimal viewing and documentation

Design meets Performance.






Mounting Options

Flexible and customized

Flexion – The name says it all, even with the different mounting options. In addition to the dimensions shown here, CJ-Optik also offers tailor-made heights and lengths of all tripod parts in order to adapt the assembly to the premises in the best possible way. In addition, the parts can be interchanged (e.g. when changing from floor to ceiling mounting).

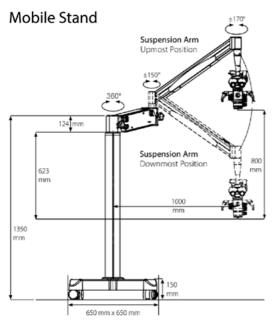


Ideal for use in multiple rooms.

Suspension Arm Upmost Position 124 mm Suspension Arm Downmost Position 1000 mm 1405 mm


The space-saving alternative to the mobile stand.

Ceiling Mount

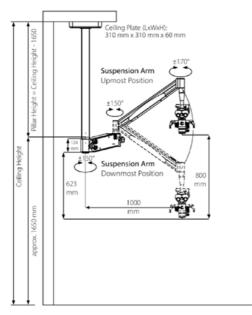

Maintains ground clearance. Pillar height can be adjusted to any ceiling height.

Wall Mount

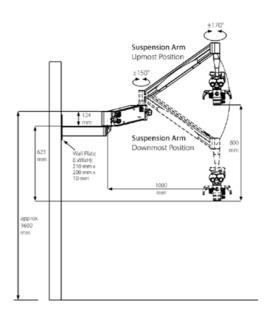

The alternative to ceiling mounting.





Optics

For **superior** results



Advanced Optical Design

CJ-Optik Flexion microscopes are highly acclaimed for their superior image quality. It is only through the communication between the optics designer and CJ-Optik that the optimal specifications can be developed and implemented. With the use of apochromatic optics in the Flexion Advanced microscope dentists can detect the finest colour and structure details.



Perfect Illumination

The fanless LED spot lighting provides a well-lit treatment area with high colour rendering index and 5.400 K. Future-proof with 65,000 hours lifespan. The spot diaphragm always shows the assistant the exact treatment site and prevents light from falling into the patient's eyes.



VarioFocus

Choose from either the VarioFocus² with 200 - 350 mm range or the VarioFocus³ with 210 - 470 mm working distance and plan-apochromatic correct built-in optics. The superior light transmission provides optimal viewing and documentation results. Removable protection cover lens with Hydrophobic Coating (HPC). Index Matching Coating for highest transmission and low reflection rates.




Details make the Difference

Just to name some of the many remarkable features that makes working with the Flexion Advanced microscopes so comfortable:

- · Individual eyecup adjustment
- Large eye relief with bigger exit pupil
- IPD adjuster optionally available

Software

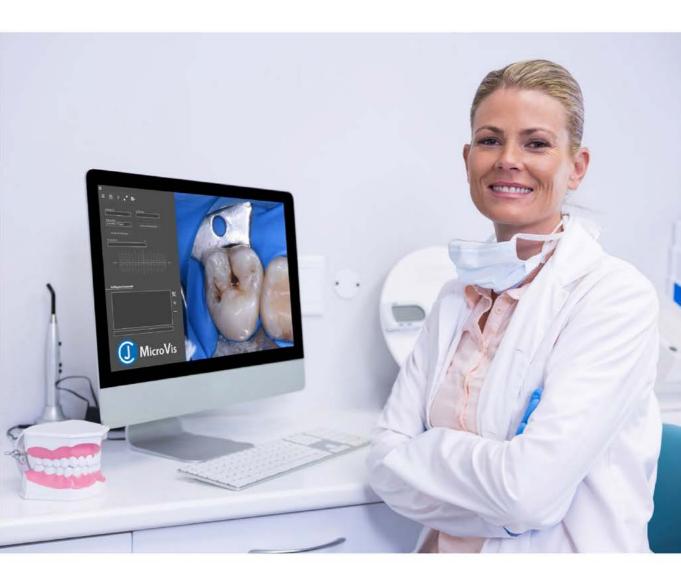

Save Time for the Important Things

With the MicroVis software specially developed for the Flexion microscopes, the image and video files of the treatment are conveniently and securely stored in your patient's file (PACS). Thanks to its user-friendliness, MicroVis relieves your team and you save important time. MicroVis is optionally available and compatible with all common practice management systems.



- Ideal for communicating with patients & assistants effectively
- · Standardizing video and image documentation workflow
- · Patient management
- · Editing of photos and videos
- Supports DICOM modalities



Documentation

Whether you prefer to document the treatment in 4K, in HD or with a smartphone: Focal lengths and optics have been uniquely designed for an ideal match to today's latest digital cameras.



4K-Imaging-PortFor Sony, Nikon & Canon
FullFrame cameras



HD-Imaging-Port For Sony, Nikon & Canon APS-C cameras



Phone-Imaging-Port

Adapter compatible
with all smartphones.
MicVision App





SensorUnit
For Sony cameras
(FullFrame & APS-C)



Full Connectivity
Swivel Arm Advanced
with HDMI, 24V, USB.
Direct power supply for the
camera from the microscope.

CJ-Optik recommends

For stress-free treatments

Ergonomics

Flexion

Anyone who has ever worked with one of our microscopes knows why we called them Flexion. We consider the Flexion as a flexible extension of the dentist's eyes. The dentist's upright treatment position ensures a relaxed posture and thus prevents neck and back problems in the long term.



MonoGlobe

The uniquely designed MonoGlobe is the weightless movement balancing system. The coupling system allows fluid and smooth repositioning at all angles without having to adjust the tension.



Quality at the Touch of a Button

Ergonomically placed controls ensure that all functions are immediately at your fingertips for all inter-procedure changes.



Cable Management

HDMI, USB, camera AC/DC power and 24V monitor power supply are integrated within the arm for a cleaner cable management. USB, HDMI, Medical monitor power integrated for a better workflow. No fragile fibre-optic cables that break or tangle.





